

IEM[®]
on life's side



CART BP_{pro}

A new era in 24-hour
blood pressure monitoring

CART BP_{pro}

Beyond the cuff

For decades, ambulatory blood pressure monitoring has been an established clinical standard around the world. Yet one thing has remained unchanged.

The cuff.

Today, a new era begins.

CART BP pro reimagines 24-hour blood pressure monitoring with a fundamentally new approach. Designed around everyday life rather than the device, it combines clinical confidence with an entirely new monitoring experience.

Not just a new device.

A new way to experience 24-hour blood pressure monitoring.




IEM[®]
on life's side

A new way to monitor blood pressure

For the first time, clinical 24-hour blood pressure monitoring is available in a lightweight ring designed for everyday life. CART BP pro uses advanced optical sensing to continuously capture physiological signals, enabling continuous 24-hour blood pressure monitoring.

The result is a redefined monitoring experience – combining clinical performance with unprecedented patient comfort.

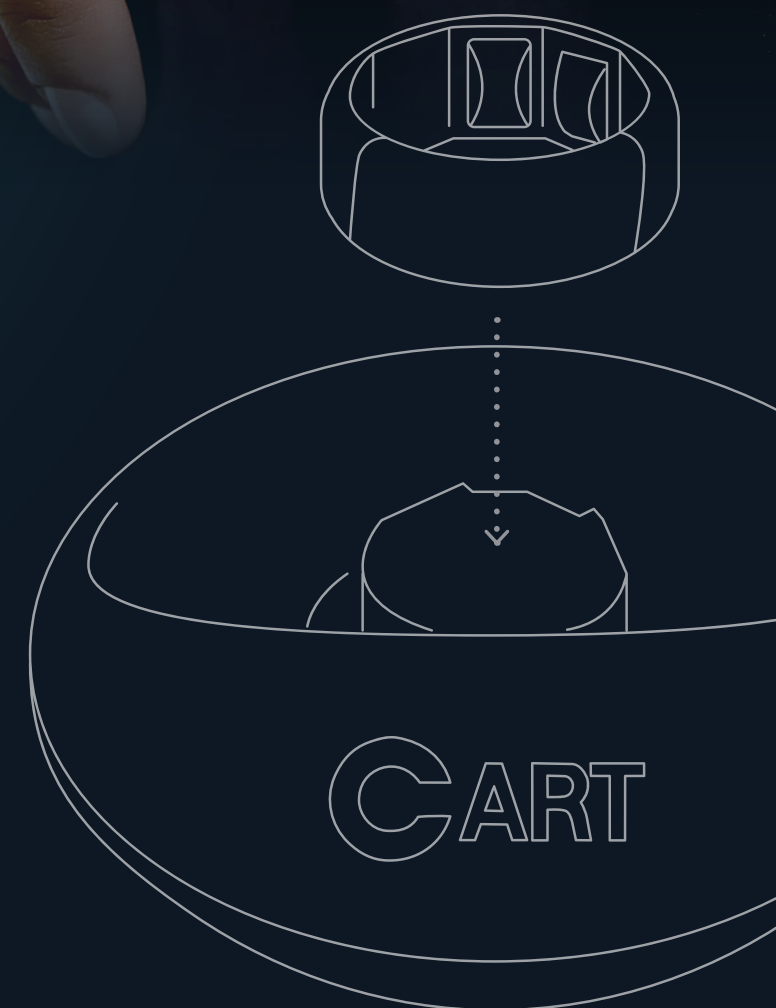



Advanced & Precise Signal Analysis

Enables precise continuous blood pressure monitoring based on physiological signals.

Photoplethysmography (PPG) Technology

Detects subtle blood flow changes in the finger vessels to accurately extract personalized blood pressure data.

Designed for Everyday Use

- ✓ Comfortable and discreet
- ✓ Convenient docking station for effortless charging
- ✓ Seamless Bluetooth® connectivity for smooth and reliable workflows
- ✓ Clinically validated measurement technology
- ✓ Continuous blood pressure insights during daily activities and sleep
- ✓ Waterproof (IP58)



↓	Patient-ID	Last Name	First Name	Date of birth	Scheduled for
	91147966	Ramirez	James	9/16/30	2/24/26, 8:00 AM
	52561665	Hanna	Elizabeth	4/2/35	3/5/26, 8:00 AM
	33565042	Flores	Joseph	6/24/81	2/28/26, 8:00 AM
	96898882	Ellison	Kimberly	1/20/74	3/2/26, 8:00 AM
	13946878	Black	Rhonda	7/6/60	3/4/26, 8:00 AM
	42188279	Reyes	Harold	3/7/35	2/24/26, 9:00 AM
	46789891	Harrison	Mary	7/23/94	3/10/26, 7:00 AM
	16817722	Douglas	Emily	10/3/50	3/4/26, 10:00 AM

↑	Type	Patient-ID	Last Name	First Name	Date of birth	Scheduled for	Serial number	Sent to device
		97819617	Andrews	Martha	10/1/66	3/1/26, 9:00 AM	753896257	3 months ago
		54138277	Moon	Howard	12/22/27	3/2/26, 7:00 AM	104494	3 months ago
		42289164	Bishop	Kelsey	1/10/28	3/2/26, 7:00 AM	753896257	3 months ago
		94447264	Day	Nicole	7/1/55	3/7/26, 7:00 AM	753896257	3 months ago
		98971620	Chavez	Linda	12/10/00	3/8/26, 7:00 AM	753896257	3 months ago

FLOW
The way to simplicity!



Validated for Clinical Confidence Evidence you can trust

CART BP pro has been clinically validated against conventional 24-hour ABPM in accordance with the validation criteria recommended by the European Society of Hypertension (ESH). The studies confirm that CART BP pro meets the defined performance criteria across 24-hour blood pressure monitoring.



Algorithm Validation with A-Line Data ¹

	Criteria [mmHg]	SBP [mmHg]	DBP [mmHg]
Mean Error	< ±5	0.98 ✓	0.52 ✓
SD	< ±8	6.97 ✓	4.38 ✓

Trained on 4,185 subjects and 320,535 A-line BP samples from an operating room, with 10% reserved for validation



Validation against Standard Sphygmomanometer ²

	Criteria [mmHg]	SBP [mmHg]	DBP [mmHg]
Mean Error	< ±5	0.16 ✓	-0.07 ✓
SD	< ±8	5.90 ✓	4.68 ✓

Demonstrated accuracy within ISO 81060-2:2018 criteria (5 ± 8mmHg) in a study of 89 participants



24-Hour ABPM Comparison: Daytime-Nighttime BP Change ³

	Criteria [mmHg]	SBP [mmHg]	DBP [mmHg]
Mean Error ± SD	< 5±8	-3.39 ± 5.86 ✓	-3.74 ± 4.72 ✓
r *	≥ 0.7	0.73 ✓	0.75 ✓

* r (The correlation coefficient) quantifies the relationship between blood pressure readings obtained from CART BP pro and the reference standard, with values ranging from -1 to 1

Streamlined and fast processes

In the often hectic environment of medical practices or clinics, every second counts. **Our software FLOW provides lightning-fast and reliable access** to patient data and device information, as well as comprehensive reports. Gone are the days of long loading times and unnecessary clicks - leaving more time for the tasks that truly matter.

Effortless usability

FLOW is a software designed for doctors and their teams to focus on their patients instead of wasting time navigating complex menus or features.

- ✓ Clear structures
- ✓ Self-explanatory workflows
- ✓ Only minimal training required
- ✓ **Ease of use at its best**

Automated, comprehensive reports

Simply connect the device to your computer to read it out and get a clearly structured report. The connected device is recognized in seconds and the PDF report is created fully automatically. **This ensures fast access to results and smooth documentation with minimal effort.**

Hassle-free integration

FLOW offers a remarkably simple and **seamless integration into existing EMR systems**, allowing you to get started immediately without any technical hurdles. After the quick installation, it operates with **unwavering stability, even under high load.** It is robust, reliable and ensures a high level of protection for sensitive data.

1. Joung J, Jung CW, Lee HC et al. Sci Rep. 2023 May 27;13(1):8605
2. Kim J, Chang SA, Park SW. J Korean Med Sci. 2024 Jan 15;39(2):e18
3. Lee H, Park S, Kwon H et al. Korean Circ J. 2024 Feb;54(2):93-104



www.iem.de

IEM[®]
on life's side

IEM GmbH
Gewerbepark Brand 42
52078 Aachen
Germany

P +49 241 41259-300
F +49 241 41259-311
customer.relations@iem.de

Linked  kununu[®]